Parasitology Stool Collection Kit

IMPORTANT: Please read and follow instructions before starting.
This collection kit will help you collect a stool sample for parasitic analysis.

1

The collection kit consists of two different vials containing liquid. DO NOT DRINK! Keep out of reach of children and pets.

2

Open kit and fill out: 1) patient name, (2) date and time specimen was collected. Please check the appropriate box for specimen description.

3

Collect the stool in a wide mouth container, bedpan, on a clean newspaper, or in a plastic bag placed over the toilet seat opening. DO NOT pass the specimen into the toilet or directly into the collection vial. Do not expose sample to urine or water.

4

Open the vials carefully. Using the collection spoon attached to the cap, add enough stool until the liquid reaches the ARROW on the label. Fill only one vial at a time and replace the cap onto the same vial. DO NOT mix caps.
IMPORTANT: Sample area of the specimen which appears bloody, slimy, or watery. If the stool is hard, sample from both ends and the middle of the specimen.


5

Ensure the cap is tightly closed and shake the tube vigorously to mix.

6

Wash hands thoroughly with warm water and soap.

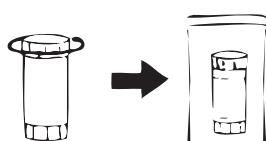7

Place the vial(s) in the original package and seal. Take specimens at room temperature to the county public health department as soon as possible.

South Carolina Department of Health
and Environmental Control

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